

Point Thumb - OI

Product Overview

www.pointdesigns.com



point designs



POINT THUMB - OI

Near MP amputation
Requires metacarpal implant

- 3 sizes: 60mm, 65mm, 70mm

Point Designs offers thumb, full finger, and partial finger mechanical prosthetic fingers that interface with osseointegrated implants. These mechanical devices offer increased range of motion and greater durability compared to standard cosmetic restorations used with osseointegrated implants.

Articulation

- Articulation at all joints
- 5 locking position in flexion

Abutment Connection

- Configurable adapter
- Clamps securely to abutment
- Adjustable rotation
- Removable via M3 screw

Usability

- Integrated durable grip covers
- Automatic spring-back at full flexion
- Lightweight

Safety


- Safety release lever
- Rotation coupler slips under high torque

Devices are intended for light-duty activities




How to Operate


Donning & Doffing

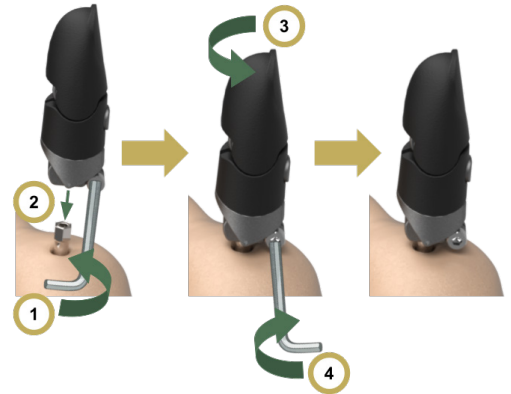

Flexion & Extension

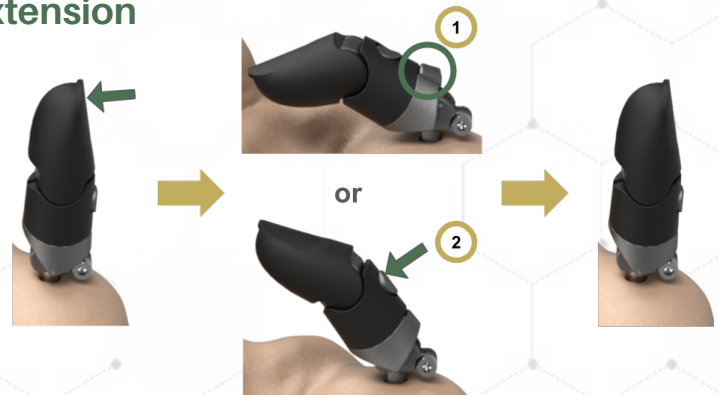

Rotation

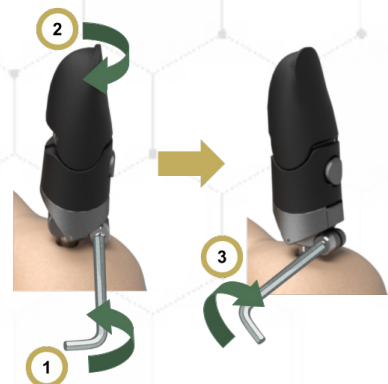

See www.pointdesigns.com/en/resources for more information